



The leader in point-of-care digital display
mounting and mobility

StyleView[®] Carts



Notebook Cart
non-powered



LCD Cart
powered

StyleView Non-Powered Carts

Lean. Clean. Smart. Cart.



1. Effortless landscape-to-portrait monitor adjustment
2. 25° of fluid independent monitor tilt
3. Side-to-side monitor pan improves data visibility for caregiver, patient and family
4. 13 cm (5") independent monitor height adjustment, for 63 cm (25") total height range, increases comfort and productivity for personalized viewing
5. 51 cm (20") of worksurface height adjustment for comfortable data entry sitting or standing with Constant Force (CF) patented motion technology
6. Lighted keyboard tray for nighttime usage
7. Ergonomic negative-tilt slide-out keyboard tray with left or right mouse tray ensures wrists remain in a neutral position during data entry—avoid repetitive-stress injuries and increase comfort and productivity
8. Electronic autolock drawer, programmable with keypad
9. Four ultra-smooth gliding casters for extremely easy push/pull mobility; two locking casters ensure lock-in-place stability



Designed for Nurses

- Take it for a spin—StyleView quickly responds to caregivers' touch and effortlessly maneuvers to the point of care
- StyleView's small footprint with ultra-smooth gliding casters travels across crowded carpeted or hard floors with ease—perfect for space-constrained areas
- Comfortable, ergonomic data entry—work sitting or standing—helps promote caregiver wellness

Engineered for IT

- Accommodates multiple configurations, including open or closed notebooks with monitor, CPU and monitor, thin clients, tablets and all-in-one computers
- Advanced open-architecture design supports future hardware upgrades
- Integrated four-port USB hub and Ethernet ready-on function: awoken enclosed computer from energy-saving hibernation mode with the touch of an external button
- Secure, ventilated storage for notebook, thin client and CPU up to 3.75" (9.5 cm) thick

Priced for Procurement

- Features include basket, front/rear handles, scanner holder, mouse pouch, advanced cable management, worksurface, four-port USB hub and Ethernet ready-on function and AC inlet powered coiled cord
- **Lowest total cost of ownership of any point-of-care cart.** StyleView is tested to UL/EN/IEC specifications through an ISO certified regimen for superior quality. Using StyleView means less maintenance and on-site technical support costs
- Extensive 5 year warranty

Certified for Safety

- We have tested to UL/EN/IEC specifications for enhanced patient safety
- Environmentally clean: RoHS, WEEE and EU (94/62/EC) packaging compliant
- Extremely easy to clean design promotes infection control

StyleView Powered Carts

Moving at maximum efficiency.



We listened! Point-of-care carts that deliver caregiver efficiency, patient safety and quality care


The power system has an unprecedented two-year warranty* which leapfrogs over warranties offered with competitive carts. Ergotron's upgraded StyleLink software enhances power system management, helping extend battery life. For patient and caregiver safety, the powered carts and power system are UL and EN 60601-1 certified. A programmable, autolock drawer facilitates medication delivery. StyleView bridges point-of-care computing and bedside barcoding, reducing healthcare costs and improving patient care.

Improve viewing comfort and enhance productivity through StyleView's expansive range of effortless height adjustment. 51 cm sit-to-stand motion is combined with an additional 13 cm independent monitor height adjustment and 25° screen tilt, which benefits all users, including those who wear progressive lenses. The keyboard tray provides up to 12° of negative tilt—recommended by ergonomists for the proper typing position. This reduces the risk of repetitive stress injuries. An enclosed storage area makes it easy to manage cables.



1. 13 cm (5"), independent monitor height adjustment increases comfort and productivity for personalized viewing, including for those wearing progressive lenses. Rotates from portrait to landscape mode. Compatible with Tablet PCs (LCD carts only)
2. Comfortable data entry, sitting or standing. Lift Range: LCD/worksurface adjust 51 cm (20") in tandem; LCD adjusts 13 cm (5") independently; maximum LCD height adjustment 63 cm (25")
3. Patented CF motion technology provides fluid, one-touch height adjustment
4. Secure, autolock drawer controlled by electronic keypad facilitates medication delivery—enhancing caregiver efficiency and patient safety. Removable cassette is configured with six bins (each with a clear lid and ID label)
5. Barcode scanner holder provides easy access to hand-held scanner
6. Wrists remain in a neutral position during data entry with negative-tilt keyboard tray—avoid repetitive-stress injuries and increase comfort and productivity
7. Expansive cable management storage area
8. Easier mousing and typing with rigid keyboard platform
9. Long-life power performance supported with an industry-leading two-year warranty* on power system. Five-year warranty on cart
10. Upgraded StyleLink™ software improves power system management tools, helping extend the life of the battery and reduce costs. This version includes email alerts for battery replacement and charge remaining. Windows Vista™ compatible (powered carts only)
11. Easy-to-clean external 3 m (10-foot) coiled cord included with powered carts
12. Easy-to-use cart comes fully assembled

* See complete warranty on website

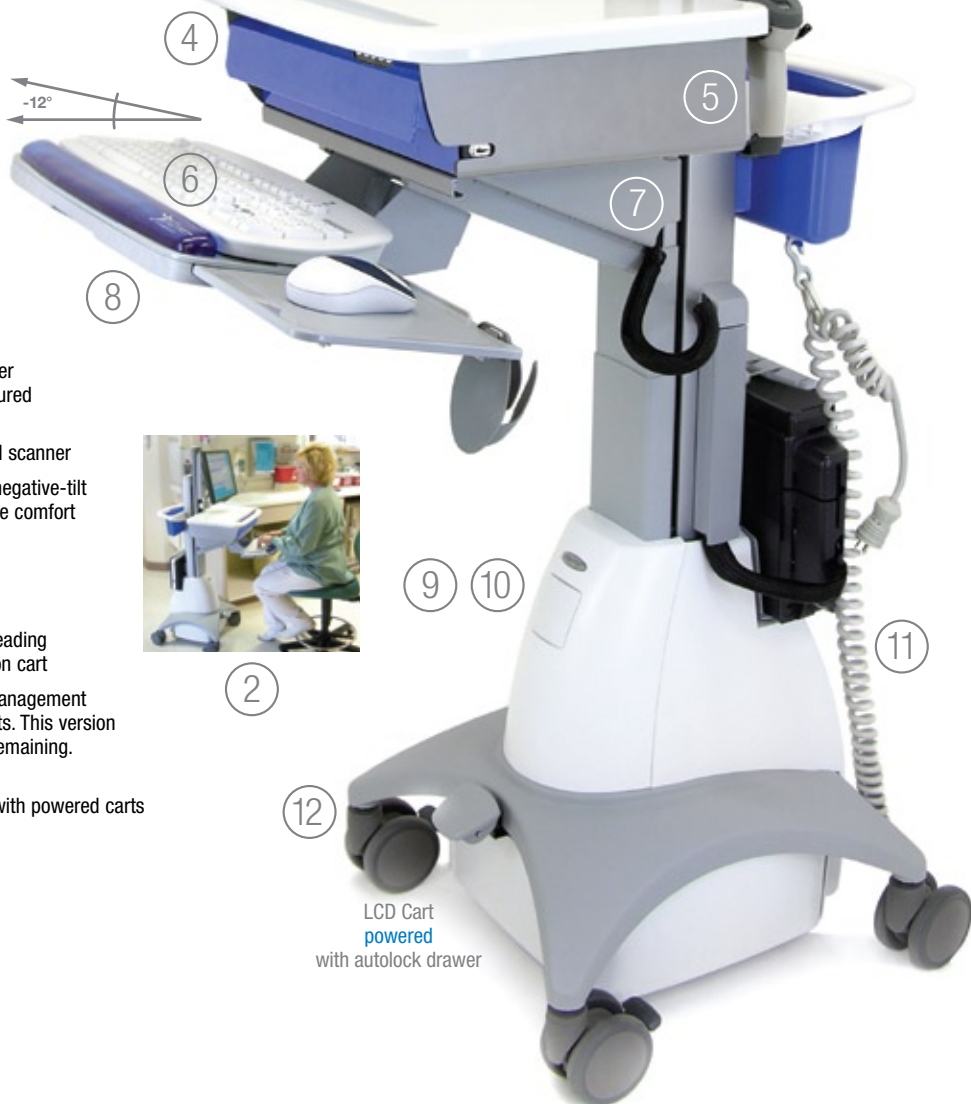

StyleView Powered Carts

Features include worksurface, four-port USB hub and Ethernet ready-on function, scanner holder bracket, mouse holder, wrist rest, negative-tilt keyboard tray, left/right mouse tray, front/rear handle, storage basket, cable storage, four casters (two locking).



LCD/Tablet PC Size	LCD Weight*	Lift **	Tilt	Pan	Rotation †	VESA
≤ 22"	2.7–7.3 kg 6–16 lbs	51 cm 20"	25°	24°	90°	MIS-D

** Lighter or heavier loads can be locked in place: < 2.7 kg (6 lbs) and 7.7–9.1 kg (17–20 lbs)

** Lift Range: LCD/worksurface adjust 51 cm (20") in tandem; LCD adjusts 13 cm (5") independently; maximum LCD height adjustment 63 cm (25")

† Portrait orientation may be incompatible with > 19" widescreen LCDs

Notebook Size*	Notebook Weight	Lift	Tilt	Pan	Rotation	VESA
≤ 17"	2.3–5.4 kg 5–12 lbs	51 cm 20"	10°	—	—	N/A

* Maximum physical dimensions of a notebook used in the cart: 43.1 x 31.7 x 3.8 cm (17" x 12.5" x 1.5")

Dual Display Size	Dual Display Weight*	Lift	Tilt	Pan	Rotation	VESA
≤ 22"	≤ 7.7 kg ≤ 17 lbs	48.3 cm 19"	30°	—	—	MIS-D

* Per LCD pivot

StyleView Non-Powered Carts

Standard cart features: worksurface (except Dual Display Cart), basket and rear handle, scanner holder, cable management system, universal CPU holder (LCD models only), full-size keyboard tray with left/right mouse tray, gel wrist rest, mouse holder, height position lock (except Dual Display Cart), high-performance casters with dual wheels for carpet or hard floor use. Powered carts include 55Ah, 12VDC, Sealed Lead Acid (SLA), Absorbed Glass Mat (AGM), Deep Cycle, Valve Regulated type battery with supporting electronics.



LCD Size	LCD Weight	Lift *	Tilt	Pan	Rotation**	VESA
≤ 22"	5–18 lbs 2.3–8.2 kg	20" 51 cm	25°	—	90°	MIS-D

* Lift Range: LCD/worksurface adjust 20" (51 cm) in tandem; LCD adjusts 5" (13 cm) independently; maximum LCD height adjustment 25" (63 cm)

** Portrait orientation may be incompatible with > 19" widescreen LCDs

Notebook Size*	Notebook Weight	Lift	Tilt	Pan	Rotation	VESA
≤ 17"	5–12 lbs 2.3–5.4 kg	20" 51 cm	10°	—	—	N/A

* Maximum physical dimensions of a notebook used in the cart: 17" x 12.5" x 1.5" (43.1 x 31.7 x 3.8 cm)

Dual Display Size	Dual Display Weight*	Lift	Tilt	Pan	Rotation	VESA
≤ 22"	≤ 17 lbs ≤ 7.7 kg	19" 48.3 cm	30°	—	—	MIS-D

* Per LCD pivot

Tablet PC Size	Tablet PC Weight	Lift *	Tilt	Pan	Rotation	VESA
≤ 22"	5–12 lbs 2.3–5.4 kg	20" 51 cm	100°	—	90°	MIS-D

* Lift Range: LCD/worksurface adjust 20" (51 cm) in tandem; LCD adjusts 5" (13 cm) independently; maximum LCD height adjustment 25" (63 cm)

	Cart	Part Number	VAC	Hz	Certification Standard	Certification Mark	Grounded Plug & Outlet
Powered							
Europe*	LCD Notebook Dual Display	SV22-92014 SV22-92011 SV22-22017	230	50	EN 60601-1	CE	European CEE 7
Ireland United Kingdom	LCD Notebook Dual Display	SV22-92015 SV22-92012 SV22-22018	230	50	EN 60601-1	CE	United Kingdom BS 1363
Hong Kong Malaysia Singapore	LCD Notebook Dual Display	SV22-92015 SV22-92012 SV22-22018	230	50	IEC 60601-1	Not Applicable	United Kingdom BS 1363
China	LCD Notebook	SV22-92024 SV22-92025	230	50	IEC 60601-1	Not Applicable	GB2009-1 GB1002-1
Australia New Zealand	LCD Notebook Dual Display	SV22-92016 SV22-92013 SV22-22019	230	50	IEC 60601-1 (AS/NZ 3200.1.0)	RCM	Australian AS/NZS 3112
Japan	LCD Notebook	SV22-92029 SV22-92030	100	50/60	JIS-T-0601	PSE	JIS 8303
Korea	LCD Notebook	SV22-92027 SV22-92028	220	60	KSC 60950	EK	European CEE 7
Saudi Arabia	LCD Notebook	SV22-92027 SV22-92028	220	60	KSC 60950	EK	European CEE 7

* Contact Ergotron Sales Representative for country-specific configuration

	Cart	Part Number
Non-powered		
Global	LCD LCD	SV31-31001
	w/autolock drawer	SV31-31051
	Notebook Notebook	SV31-31002
	w/autolock drawer	SV31-31052
	Dual Display	SV21-31645



© 2008 Ergotron, Inc. rev. 21/05/2008 **Content is subject to change without notification**
Some products may not be available in some regions of the globe—visit our website for more information

Americas Sales and Corporate Headquarters	EMEA Sales	APAC Sales	Worldwide OEM Sales
St. Paul, MN USA (800) 888-8458 +1-651-681-7600 www.ergotron.com sales@ergotron.com	Amersfoort, The Netherlands +31 33 45 45 600 www.ergotron.com info.eu@ergotron.com	Singapore www.ergotron.com info.apac@ergotron.com	www.ergotron.com info.oem@ergotron.com

