

# ergotron<sup>®</sup>

*moving you forward*

- Helps reduce cart downtime with real-time insights into usage patterns, drawer access, battery life and more through an online dashboard
- Software communicates directly with the cloud without access to patient records to keep information secure and support HIPAA compliance
- Provides real-time notifications to your IT staff when a mobile workstation requires service or attention to eliminate guesswork and unnecessary checks
- Enables IT staff to diagnose issues from a distance to help reduce support time instead of traveling across expansive campuses to locate a cart
- Based on your service level agreement, automated service notifications allow Ergotron to address an event as it occurs to remotely diagnose the situation and dispatch service and parts
- Free with your StyleView 42 or StyleView 44 medical cart and quick to implement with installation per device in as little as one minute

eKinnex<sup>™</sup>  
SOFTWARE

FOR  
STYLEVIEW<sup>™</sup>

## CLOUD-BASED FLEET MANAGEMENT

21st century, patient-centered healthcare needs 21st century technology and equipment. Rethink how you manage your medical cart fleet with eKinnex<sup>™</sup> technology. This secure, web-based software hosted by Ergotron monitors your fleet in real time, helping you proactively identify service and maintenance needs to maintain uptime, freeing up time for your IT team. The HIPPA-compliant portal keeps patient information safe while providing key insights to manage your fleet to meet the needs of fast-paced, modern care.

## Flexible access control and cart configuration

Control who can access drawers through the cart users page. For more advanced access, separate carts into groups and control access on a group basis, along with settings like mute drawer alarm and open drawers at once.

## Detailed power system metrics



Historical power system data, including charge, load and current, lets you quickly identify usage patterns and potential issues.

## Drawer access log

Use the drawer log page to view and search drawer events to see who opened which drawer on which cart at what time.

## Easy to deploy

The cart client software installs on your cart-connected computers using a Microsoft Installer package and minimal resources—just 2 MB of memory. To get started, follow the link in your welcome email, download the software onto your local client PC and run the installation file.

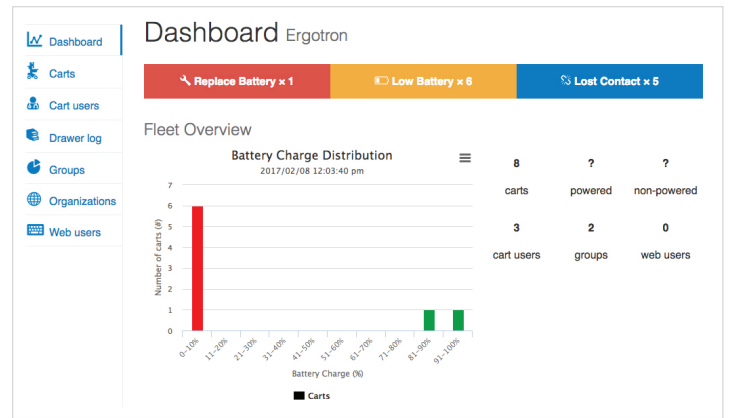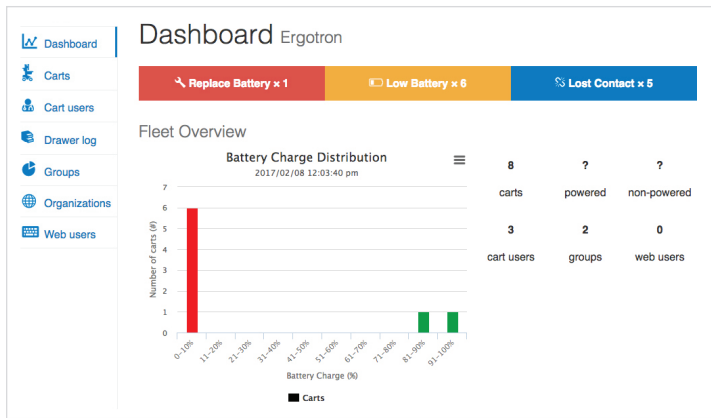

**Carts Ergotron**

Replace Battery × 1 | Low Battery × 6 | Lost Contact × 5

**Filters**

Search: Name or Group

Connection:  Online,  Offline

Power Source:  Plugged In,  Battery,  None

Battery Condition:  Good,  Replace

Charge: 0% to 100% (< 100%)

**Carts (8)**

Name	Last Seen	Group	Type	Power Source	Charge	Battery Condition
ERG-TANGENT-PC	2017/01/27 3:30:18 pm	HIMSS17	Unknown	Unknown	0%	Checking
ERGLOANER08	2017/01/24 4:20:51 pm	demo	Unknown	Valence 12RTL 40AH	0%	Checking
NHAZZARD1	2017/02/08 12:06:40 pm	demo	Unknown	Unknown	0%	Checking
SOLSON-TEST	2017/02/08 12:05:57 pm	demo	Unknown	Unknown	81%	Checking
SSQUIRES1	2017/02/02 12:35:44 pm		Unknown	Unknown	0%	Checking
SURFACEPRG3	2017/02/08 8:51:57 am	demo	Unknown	Unknown	100%	Replace
ITOWNES	2017/02/02 4:36:58 pm		Unknown	Unknown	0%	Checking
USER-PC	2017/01/27 3:30:31 pm	HIMSS17	Unknown	Unknown	0%	Checking

Key features:

- Create and delete users
- Manage user information and drawer PINs
- Manage user groups
- Sort and filter to identify maintenance needs
- Import and export users

Learn more at [www.ergotron.com/ekinnex](http://www.ergotron.com/ekinnex).

At-a-glance insights into your fleet:

- Real-time notifications when a mobile workstation requires service or attention
- Battery status
- Last report to the cloud for the cart, power system and drawers
- Power and drawer event log
- Cart software and hardware information, like last reported location

Download additional resources at [ergotron.com](http://ergotron.com).

For more information:

Americas: 800.888.8458 / +1.651.681.7600 / [insidesales@ergotron.com](mailto:insidesales@ergotron.com)

EMEA: +31.33.45.45.600 / [info.eu@ergotron.com](mailto:info.eu@ergotron.com)

APAC: [apaccustomerservice@ergotron.com](mailto:apaccustomerservice@ergotron.com)

Custom: [custom@ergotron.com](mailto:custom@ergotron.com)

© Ergotron, Inc. 03.31.2020

Literature made in the U.S.A. Content subject to change.

Ergotron devices are not intended to cure, treat, mitigate or prevent any disease.

Patent information available at [www.ergotron.com/patents](http://www.ergotron.com/patents).

WorkFit, StyleView, LearnFit, Neo-Flex, PowerShuttle, LiFeKinnex, CareFit

and eKinnex are trademarks of Ergotron.

